

Albion products bear the **CE** mark denoting their compliance with the safety requirements of the European MDD 93/42/EEC.

Albion carbon and stainless surgical blades comply with the standards BS 2982: 1992, ISO 7740: 1985 and EN 27740:1992. Our range of blade sizes are 10, 10A, 11, 12, 15, 20, 21, 22, 23 & 24. Blade are packed individually in an aluminum foil packet. Carbon steel blades are VPI protected. Unit containers consist of 100 blades in a colour coded dispenser carton.

The Albion Disposable Scalpel consists of a Stainless Steel blade fitted on an ABS plastic handle. The blade edge is protected by a cap to prevent pre or post surgery injury. Disposable Scalpels are individually packed in a see through peel-pouch. Colour coded unit cartons contain 10 scalpels.

All products are sterilized by Gamma Radiation to a minimum dosage of 25 kGy (2.5 Mrads).




Albion Produkte tragen das **CE** Zeichen, das bezeichnet, daß sie die europäischen MDD/93/42/EEC Sicherheitsbestimmungen einhalten.


Albion chirurgische Klingen aus Kohlenstoff- und rostfreiem Edelstahl halten die Normen BS 2982: 1992, und ISO 7740:1985 und EN 27740:1992 ein. Unsere Reihe der Klingen umfaßt die Größen 10, 10A, 11, 12, 15, 20, 21, 22, 23 und 24. Die Klingen werden einzeln in Päckchen aus Aluminiumfolie verpackt. Klingen aus Kohlenstoffstahl sind VPI geschützt. Ein Container besteht jeweils aus 100 Klingen in einem farblich kodierten Spender.

Das Albion Einmalskalpell besteht aus einer rostfreien Edelstahl Klinge, die an einem ABS Plastikgriff befestigt ist. Die Klingenschneide hat eine Schutzkappe, die eine vor- oder nachchirurgische Verletzung verhindert. Skalpelle werden einzeln in durchsichtigen Schälbeuteln verpackt. Ein farblich kodierter Karton enthält 10 Skalpelle.



Die Produkte werden mittels einer Mindestdosis von 25 kGy (2.5 Mrads) Gammastrahlung sterilisiert.

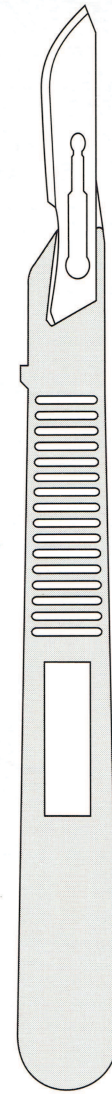

FITMENT NO: 3

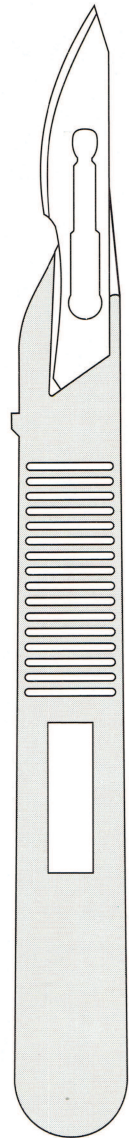

FITMENT NO: 4

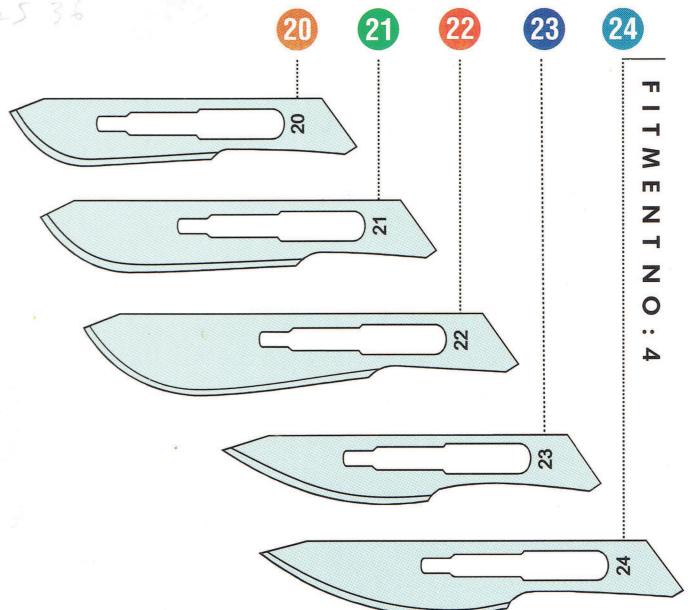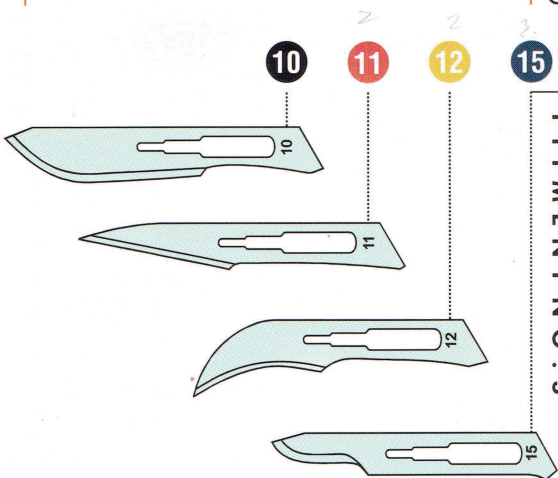

**ALBION**

Albion Surgicals Limited, Aizlewood's Mill, Nursery Street, Sheffield S3 8GG, U.K.  
Tel: (0)114 282 3310 Fax: (0)114 282 3150 e-mail: 106070.1205@compuserve.com