


# PROBEL OMEGA IV

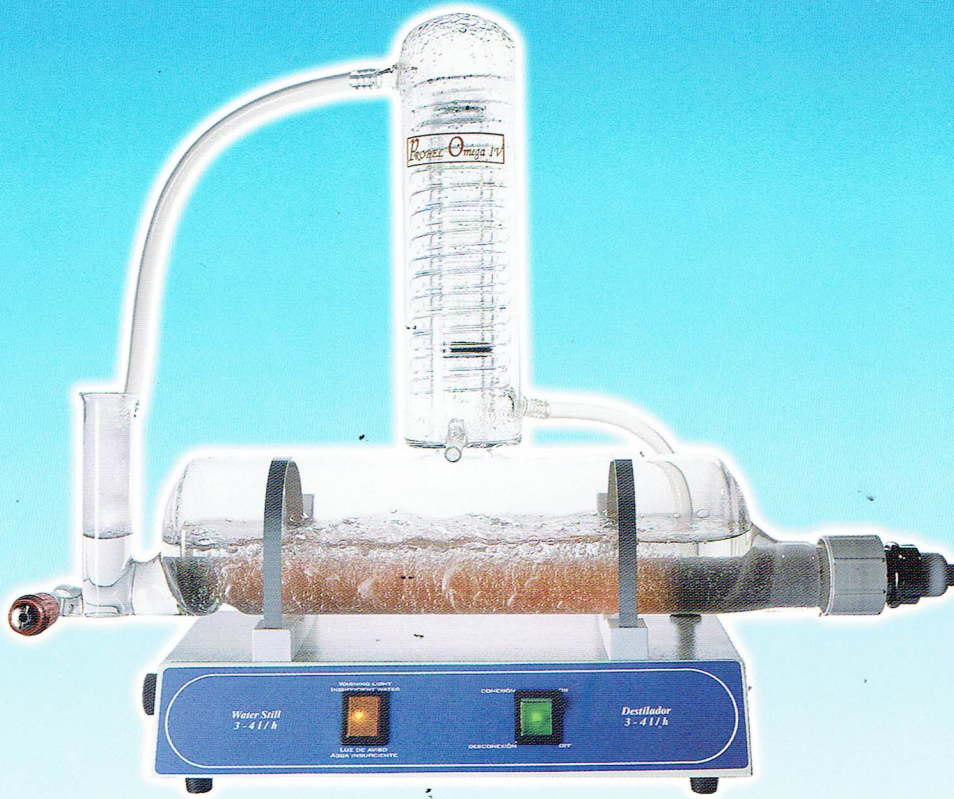

## *PROBEL Omega IV*

is designed to provide reliable continuous operation. The user friendly features are one of the many, unmatched by many. It is easy to install, easy to use, easy to maintain with safety features, is all you require of a Laboratory Water Distiller.

### FEATURES

Capable of producing 3 to 4 liters of pyrogen free distilled water.

### CONSTRUCTION

Epoxy coated stand with pressure switch and indicator lamps. Borosilicate glass condenser and Boiling chamber fitted with Silica sheathed heating element.

### SPECIFICATION

**Dimension:** 600 x 180 x 480mm  
(L x B x H)

**Weight:** 7kgs gross / 5kgs Nett

**Power:** 2,200 watts (10Amp)

**Electrical Compliance:** ESCHLE

**Water Supply:** Approx. 1 lit/min (5 P.S.I.)

### DISTILLATE QUALITY

**Resistivity:** 0.5 – 1.0 mΩ/cm

**Conductivity:** less than 2.0 μs/cm

**pH Level:** 5.0 to 6.5

**Pyrogen free:** Yes (however regular testing of your output is advised)

*Your preferred distributor*